



Proven technology. More individuality.

# Elmasonic Select

**ДИА•М**  
современная лаборатория

[www.dia-m.ru](http://www.dia-m.ru)  
заказ on-line



## Areas of use

**Powerful and individual - these are the ultrasonic cleaning devices of the Elmasonic Select series.**

**Equipped with 37 kHz ultrasonic frequency, 5 ultrasonic modes, the possibility to save up to 4 cleaning programmes, a heating and tank volumes from 3 to 90 litres (from Select 60 with drain), Elmasonic Select devices are versatile.**

**Application-specific cleaning chemicals and practical accessories such as baskets, noise protection box, holders, etc. ensure the best cleaning results and make cleaning easier.**

### Developed for

Research and development laboratories •  
Pharmaceutical and industrial-analytical laboratories • Work and repair shops • Metal processing companies • Manufacturers of industrial components • Watchmakers • Jewellers • Watch and jewellery workshops



### Cleans

Laboratory instruments such as test tubes, glass pipettes and burettes • Industrial components and parts such as carburettors, filters • Watch components such as hinges, clasps, housings • Jewellery such as earrings, bracelets - with and without gem stones



### Removes

Resins • Oil and grease • Adhesive residues • Grinding and polishing pastes • Chips • Flux agents • Organic contamination such as horn and skin residues • Wear dirt

### Prepares samples

Degassing • Dispersing • Emulsifying • Homogenising • Deagglomerating • Thermosensitive incubations • Lyses • Biochemical assays





## Individual



### Adjustable cleaning programmes

- Memory for 4 individual cleaning programmes
- Cleaning process in the ultrasonic cleaning device can be called up for reproducible results
  - Fast switching between favourite programmes



### To configure

The auto-degas mode removes disturbing air bubbles in the cleaning liquid and thus leads to better cleaning results. For temperature-sensitive parts or temperature-sensitive processes, the limit temperature can be set individually and the display changed to °C / °F.

## Strong



### 5 ultrasonic modes for best cleaning results and easy sample preparation

- dynamic:** Combines the ultrasonic modes pulse and sweep to increase the overall cleaning performance
- pulse:** Removes even stubborn contaminations through a pulsating sound field distribution
- sweep:** Uniform distribution of the ultrasonic power for most cleaning tasks

- eco:** For gentle cleaning with quiet operation and for a longer service life of the ultrasonic cleaning device
- degas:** Quick degassing of freshly prepared cleaning solutions and (HPLC) samples

## Intuitive



### Clear and simple

Clear, easy-to-read display showing all relevant parameters e.g. remaining cleaning time, set/actual temperature for precise adjustment.



### With feedback

The user is always informed about the operating status of the device via clear symbols for e.g. end of cleaning or overtemperature as well as acoustic and visual feedback.

## Application oriented




### With matching accessories

Plastic lids ensure quieter operation and serve as drip trays for baskets and parts. Further accessories such as holders, inserts or noise protection boxes facilitate the work for the most diverse cleaning tasks.



### Through industry-specific cleaning chemistry

The industry-specific cleaning chemistry is ecologically and economically optimised. Overviews for the selection of suitable cleaners with application and dosing aids are available.

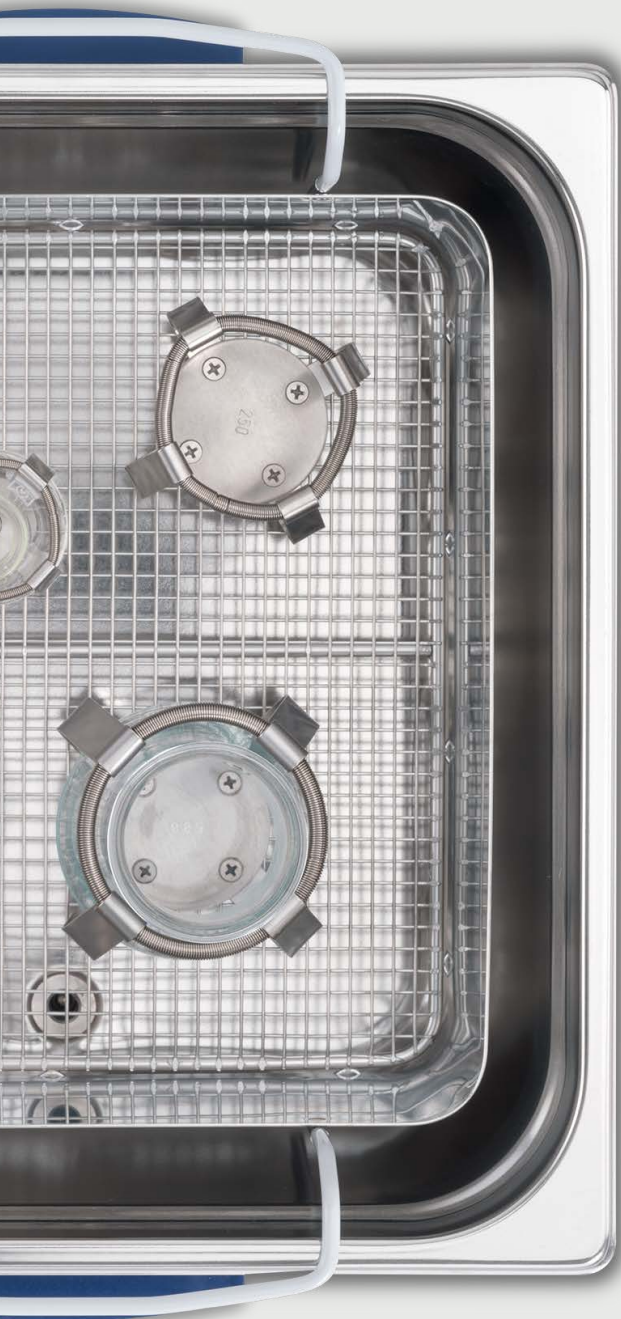

Elma Schmidbauer GmbH

**000 «Диаэм»**

Москва  
ул. Магаданская, д. 7, к. 3 ■ тел./факс: (495) 745-0508 ■ sales@dia-m.ru

[www.dia-m.ru](http://www.dia-m.ru)

**С.-Петербург**  
+7 (812) 372-6040  
spb@dia-m.ru

**Новосибирск**  
+7(383) 328-0048  
nsk@dia-m.ru

**Воронеж**  
+7 (473) 232-4412  
vrn@dia-m.ru

**Йошкар-Ола**  
+7 (927) 880-3676  
nba@dia-m.ru

**Красноярск**  
+7(923) 303-0152  
krsk@dia-m.ru

**Казань**  
+7(843) 210-2080  
kazan@dia-m.ru

**Ростов-на-Дону**  
+7 (863) 303-5500  
rnd@dia-m.ru

**Екатеринбург**  
+7 (912) 658-7606  
ekb@dia-m.ru

**Кемерово**  
+7 (923) 158-6753  
kemerovo@dia-m.ru

**Армения**  
+7 (094) 01-0173  
armenia@dia-m.ru

